

Goodman HTS Pipeline User Manual

version 0.1

Simón Torres

August 08, 2017

Contents

Introduction	1
Overview	1
Features	1
Ways to run the pipeline	1
What the pipeline does not do	2
General Considerations on using the pipeline	3
Command line arguments	3
Lists of Reference Lamps Available	4
Adding new reference lamps	5
Running the pipeline in the SOAR data reduction computer	5
Establish a VNC connection	5
VNC from the Terminal	6
VNC using a Graphical Client	6
Running the Pipeline	6
Using the Interactive Mode	8
Troubleshooting	15
Installation Instructions	16
Install Dependencies	16
Ubuntu 16.04	16
CentOS 7	16
Installing on MacOSX	16
Install Using Virtual Environments	16
Downloading the Goodman HTS Spectroscopic Pipeline	17
Installing the Pipeline	17
Install DCR	17
Install binary DCR	17

Introduction

This is the User Manual for the *Goodman Spectroscopic Data Reduction Pipeline*. It provides an overview of the pipeline's main features, instructions on its use and how to run it on our dedicated *SOAR Data Reduction Server*, and installation instructions for those who wish to run it on their own computers.

Overview

The Goodman Spectroscopic Data Reduction Pipeline - GOODSPEC - is a Python-based package for producing science-ready, wavelength-calibrated, 1-D spectra. The goal of **goodspec** is to provide SOAR users with an easy to use, very well documented software for reducing spectra obtained with the Goodman spectrograph. Though the current implementation assumes offline data reduction, our aim is to provide the capability to run it in real time, so 1-D wavelength calibrated spectra can be produced shortly after the shutter closes.

The pipeline is primarily intended to be run on a data reduction dedicated computer. Instructions for running the software are provided in the [Using Pipeline](#) section of this guide. The Goodman Spectroscopic Data Reduction Pipeline project is hosted at GitHub at [it's GitHub Repository](#).

Currently the pipeline is separated into two main components. The initial processing is done by `redccd`, which trims the images, and carries out bias and flat corrections. The spectroscopic processing is done by `redspec` and carries out the following steps:

- Identifies multiple targets (spectra of more than one object in the slit)
- Trace the spectra
- Extract the spectra
- Estimate and subtract background
- Find the wavelength solution. Defaults to automatic wavelength solution, but can be done interactively
- Linearize data (resample)
- Write wavelength solution to FITS header
- Create a new file for the wavelength calibrated 1D spectrum

Features

- Self-contained, full data reduction package for the most commonly used spectroscopic setups with Goodman. Given the almost limitless number of possible configurations available with the Goodman instrument, only the most popular configurations will be supported, though we will try to add as many modes as possible.
- Python based, using existing Astropy libraries as much as feasible.
- Extensively documented, using general coding standards: PEP8 – Style Guide, PEP257 – Docstrings Convention (in-code documentation) – Google Style
- Multiplatform compatibility (tested on Linux Ubuntu, CentOS and MacOSX).
- Modular design. Could be used as a library within other Python applications.

Ways to run the pipeline

There are two ways to use the pipeline.

1. **Run it directly on a SOAR data reduction server** that you can access using VNC.
2. **Download and install the pipeline** (go to the [Install](#) section of this manual). Though we will try our best to provide answers to quick and simple installation issues, we cannot provide general installation support.

What the pipeline does not do

- In its current version the pipeline does not perform combination of individual spectra. If you obtained several individual exposures of the same object, they will be output as separate 1-D, wavelength-calibrated spectra
- There is yet no flux calibration. We are working on a module that will do this.
- This pipeline does not evaluate nor select data by quality. It will simply try to run using all existing files. **Make sure you only have good data in the folder that will be reduced.**

General Considerations on using the pipeline

The Goodman Spectroscopic Pipeline is meant to work as a single package. However, the full process is split in two separate modules: `redccd` and `redspec`. The first does the basic 2D image reduction, applying bias and flat field corrections, and cosmic ray removal. The second module, `redspec`, takes the corrected 2D images output by `redccd` and produces wavelength-calibrated 1D spectra.

The pipeline is run from the command line in a terminal window. Each module is run separately, first `redccd` followed by `redspec`, however, you could run both sequentially from e.g. a shell script.

In order to make things easier you should organize your data:

1. Make sure all the data in your folder corresponds to the same binning, readout mode, region of interest (ROI), and grating/wavelength mode combination.
2. You should have bias, flats (quartz or dome flats), and the appropriate comparison lamps. Other files like acquisition images, slit images and focus images should be deleted.
3. Do not mix dome flats with quartz lamp flats. As an example, suppose I took both quartz lamps and dome flats for my targets. I could create two folders, one with the science data and the dome flats, and another with the same science data and the quartz lamps. Then, if I run the pipeline in each folder I can compare the results and decide which type of flat works best for my particular case.

Command line arguments

For a list of the options and command line arguments type `--help` argument:

For `redccd`

```
usage: redccd [-h] [--auto-clean] [--cosmic <method>]
               [--dcr-par-dir <dcr.par_directory>] [--debug]
               [--flat-normalize <normalization_method>]
               [--flat-norm-order <order>] [--ignore-bias] [--ignore-flats]
               [--keep-cosmic-files] [--log-file <log_file>]
               [--raw-path <raw_path>] [--red-path <red_path>]
               [--saturation <value>]

Goodman CCD Reduction - CCD reductions for Goodman spectroscopic data.

optional arguments:
-h, --help            show this help message and exit
--auto-clean          Automatically clean reduced data directory
--cosmic <method>    Clean cosmic rays from all data. Options are: 'dcr',
                     'lacosmic' or 'none'. Default is 'dcr'. See manual for
                     full description of dcr.
--dcr-par-dir <dcr.par_directory>
                   Directory of default dcr.par file
--debug              Show detailed information of the process.
--flat-normalize <normalization_method>
                   Choose a method to normalize the master flat
                   forspectroscopy. Choices are: mean, simple (model) and
                   full (fits model to each line).
--flat-norm-order <order>
                   Defines the order of the model to be fitted. Default
                   to 15
--ignore-bias         Ignore bias correction
--ignore-flats        Ignore flat field correction
--keep-cosmic-files   After cleaning cosmic rays with dcr, do not remove the
                     input file and the cosmic rays file.
--log-file <log_file>
                   Name for log file. Default name is goodman_ccd.log.
                   The file is written in <red_path> and will be deleted
                   each time you run this program
```

Lists of Reference Lamps Available

```
--raw-path <raw_path>
          Path to raw data.
--red-path <red_path>
          Path to reduced data.
--saturation <value> Saturation limit. Default to 65.000 ADU (counts)
```

And for `redspec`

```
usage: redspec [-h] [--data-path <Source Path>]
                [--proc-path <Destination Path>]
                [--search-pattern <Search Pattern>]
                [--output-prefix <Out Prefix>] [--extraction <Extraction Type>]
                [--reference-files <Reference Dir>] [--interactive] [--debug]
                [--log-to-file] [--max-targets <max targets>] [--save-plots]
                [--plot-results]
```

Extracts goodman spectra and does wavelength calibration.

optional arguments:

```
-h, --help            show this help message and exit
--data-path <Source Path>
          Path for location of raw data. Default <./>
--proc-path <Destination Path>
          Path for destination of processed data. Default <./>
--search-pattern <Search Pattern>
          Pattern for matching the goodman's reduced data.
--output-prefix <Out Prefix>
          Prefix to add to calibrated spectrum.
--extraction <Extraction Type>
          Choose a which extraction to perform. Simple is a sum
          across the spatial direction after the background has
          been removed. Optimal is a more advanced method that
          considers weights and profilefitting.
--reference-files <Reference Dir>
          Directory of Reference files location
--interactive         Interactive wavelength solution. Disabled by default.
--debug              Debugging Mode
--log-to-file        Write log to a file
--max-targets <max targets>
          Maximum number of targets to be found in a single
          image. Default 3
--save-plots         Save all plots in a directory
--plot-results       Show wavelength calibrated spectrum at the end.
```

Lists of Reference Lamps Available

The automatic wavelength calibration relies on having previously calibrated reference lamps obtained in the same configuration or mode. It is also important that the lamp names are correct, for instance `HgAr` is quite different than `HgArNe`. For interactive wavelength calibration, reference lamps are used as a visual aid only. It lets you find the matching laboratory lines values that will be used to fit a pixel to wavelength relation that we call *Wavelength Solution*. The list of lamps is the following.

Grating	Mode	Filter	Lamp
400	M1	None	HgAr
400	M1	None	HgArNe
400	M2	GG455	HgAr
400	M2	GG455	HgArNe

600-old	Blue	None	HgAr
600-old	Blue	None	CuHeAr
1200	M2	None	CuHeAr
1200	M3	None	CuHeAr
1200	M5	GG455	CuHeAr

Adding new reference lamps

It is possible to add new lamps very easily you just need a raw lamp that meets the following specifications with respect to your science project:

- Same instrument configuration or mode
- Same Grating
- Same order blocking filter if present
- Same binning
- Same lamp/combination that you use in your observations
- Smallest slit possible. Equal is OK too.

Then you can use the interactive mode or other software (such as IRAF) to produce a wavelength-calibrated 1D spectrum. Now you have to options, identify the system folder where the lamps that come with the package are saved and simply put it there or put it in another directory and use the argument --reference-files

```
redspec --reference-files /path/to/ref-lamp-location
```

Or send it to me and I will make it available as a package filea.

Running the pipeline in the SOAR data reduction computer

The Goodman Spectroscopic Data Reduction Pipeline has been installed on a dedicated computer at SOAR. The procedure is to open a VNC session, for which you need to be connected to the SOAR VPN. The credentials for the VPN are the same you used for your observing run, provided by your *Support Scientist*, who will also give you the information for the data reduction computer VNC connection.

Establish a VNC connection

Separately, you should receive a server hostname, IP, port and VNC-password. If you don't you can ask for it. We have decided to use a similar organization of vnc desktops:

VNC ports and working folder assigned to each partner.

Port	Partner/Institution	Folder
:1	NOAO	/home/goodman/data/NOAO
:2	Brazil	/home/goodman/data/BRAZIL
:3	UNC	/home/goodman/data/UNC
:4	MSU	/home/goodman/data/MSU
:5	Chile	/home/goodman/data/CHILE

For the rest of this tutorial we will assume your host name is `vnc-server` the port is 1 and your password is `password`. Though we recommend using RealVNC, most other VNC clients will work fine (e.g., Remmina in Linux). For GNU/Linux and Mac OSX machines we suggest the RealVNC Viewer client. For Windows machines, we suggest

Running the Pipeline

either the RealVNC Viewer client or the UltraVNC viewer client. We also know that Vinagre and vncviewer on GNU/Linux work fine.

VNC from the Terminal



Open a terminal, and assuming you have installed vncviewer.

```
vncviewer vnc-server:1
```

You will be asked to type in the *password* provided.

VNC using a Graphical Client

Using a graphical VNC client is quite similar and intuitive

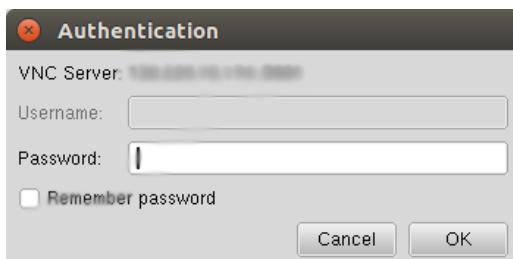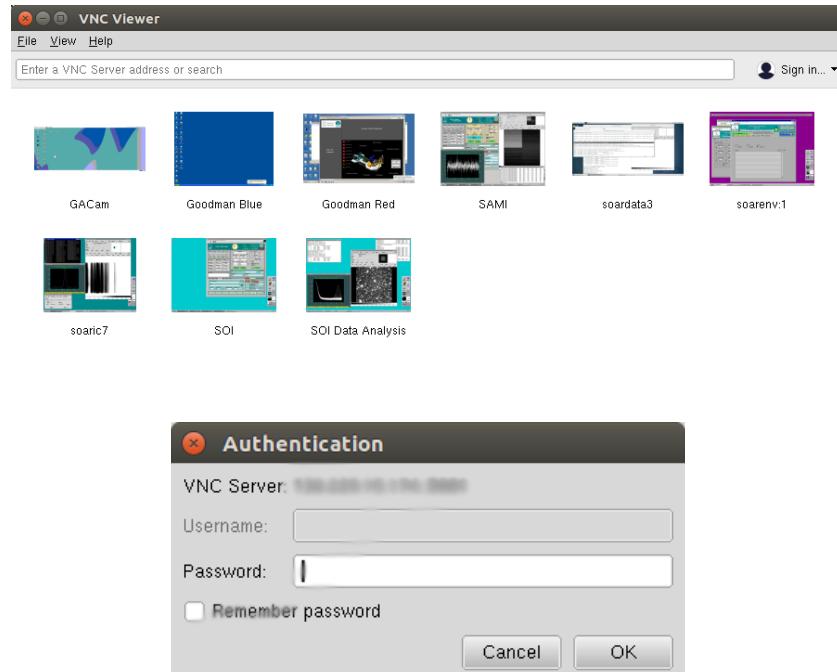

In this case the *IP address* was used, which is equivalent and sometimes better.

Running the Pipeline

1. Open a Terminal

2. Go to /home/goodman/data

```
cd /home/goodman/data
```

3. Here you have a workspace to put your data according to your institution.

The image shows a terminal window with a dark background. The window title is 'goodman@soardata3:~/data'. The terminal prompt is '[goodman@soardata3 ~]\$'. The user runs several commands: 'pwd' (prints the current working directory), 'cd data' (changes the directory to 'data'), 'ls' (lists files in the 'data' directory), and 'df -h' (displays disk usage). The output of 'ls' shows files named BRAZIL, CHILE, MSU, NAO, test, and UNC. The output of 'df -h' is as follows:

Filesystem	Size	Used	Avail	Use%	Mounted on
/dev/nvme0n1p1	24G	11G	12G	49%	/
devtmpfs	16G	0	16G	0%	/dev
tmpfs	16G	172K	16G	1%	/dev/shm
tmpfs	16G	18M	16G	1%	/run
tmpfs	16G	0	16G	0%	/sys/fs/cgroup
/dev/nvme0n1p2	211G	61M	201G	1%	/home_local
/dev/md127	5.5T	42G	5.4T	1%	/home
tmpfs	3.2G	20K	3.2G	1%	/run/user/8143
tmpfs	3.2G	0	3.2G	0%	/run/user/8142

Running the Pipeline

4. Create a data folder inside your workspace.

```
cd NOAO  
mkdir 2017-07-05  
cd 2017-07-05
```

5. Copy your data from Goodman Computer

```
scp observer@soaric7:/home3/observer/GODMAN_DATA/NOAO/2017-07-05/ ./
```

6. Make sure you have a full data set. At this point your observing logs will become very useful, eliminate focus sequence, aquisition exposure and any other file present that will not be needed for the processing. The following list summarizes the kind of data that you need to fully process your data.

- BIAS: Bias
- FLAT: Flats
- COMP: Comparison Lamps
- OBJECT: Science Frames

Also make sure your data has the same *readout speed*, *binning*, and *ROI*. If you used different configurations during the same night, we recommend you to set up a separate folder for each.

7. Run redccd:

For `redccd` I suggest using `--cosmic` and `auto-clean` also you might want to consider `--saturation <new value>` to change the saturation level if you get all your flats rejected due to saturation. Sometimes there is a hot column at the end that produced very high values.

```
redccd --cosmic --auto-clean
```

In case you want to use `--saturation` here is an example:

```
redccd --cosmic --auto-clean --saturation 70000
```

This changes the saturation level to *70000 ADU* in this context the saturation value works as a threshold for rejecting images.

By default, `redccd` puts reduced data in a subdirectory `RED`, you can provide a different one by using `--red-path`.

An image `image_file.fits` that has been fully (and propperly) processed should have the new name (including the reduced data folder):

```
cfzsto_image_file.fits
```

The meaning of every letter in `cfzsto` is summarized in the following table:

Letter	Meaning
c	Cosmic ray cleaned or mask created depending on the method
f	Flat corrected
z	Zero or Bias corrected
s	Slit trimmed, trims off the non-illuminated sections of the detector
t	Initial Image trimming
o	Overscan corrected

8. Run redspec:

By default `redspec` will search for images with the prefix `cfzsto`, in case you have produced a different prefix you can change it by using `--search-pattern`

You can just run `redspec` in case everything is the default but if this is the first time you run the pipeline I suggest:

```
redspec --plot-results
```

In that way two important plots will be shown full screen, the comparison lamp fitted to a reference comparison lamp and some values for the wavelength solution fit and the extracted spectrum plotted with the wavelength solution.

The final image has a `g` added to the start of the name, following the above example your final 1D and wavelength calibrated image will be named:

```
gcfzsto_image_file.fits
```

Using the Interactive Mode

If you need to make sure that your solution is the best possible you can use the *Interactive Mode*. Some of the key features are:


- Manually match spectroscopic lines
- Zoom in to get a better reference
- Evaluate the quality of your solution
- Uses a visual reference for matching the lines
- Uses laboratory values for reference lines

Right now it uses matplotlib as the underlying tool to enable interaction.

Below you will see a sequence as well as a description of the procedure to use the interactive mode.

1. The interactive screen has four subplots. In the *Upper main panel* you have the raw comparison lamp, i.e. your comparison lamp 1-D-extracted but without wavelength solution, the dashed red lines represent the lines found by the pipeline. In the *Lower main panel* you have the reference lamp. This is a previously calibrated lamp that is distributed with the package (also you can use your own). In this case the red lines are the laboratory values obtained from the NIST Atomic Spectra Database site (https://physics.nist.gov/PhysRefData/ASD/lines_form.html). The *Upper right panel* is a static help, intended to give you a quick and easy-to-reach help. Finally, the *Lower right panel* is an information window, in fact it serves three main purposes:

- Display messages and warnings.
- Shows you a zoomed line.
- Displays the wavelength solution quality information.

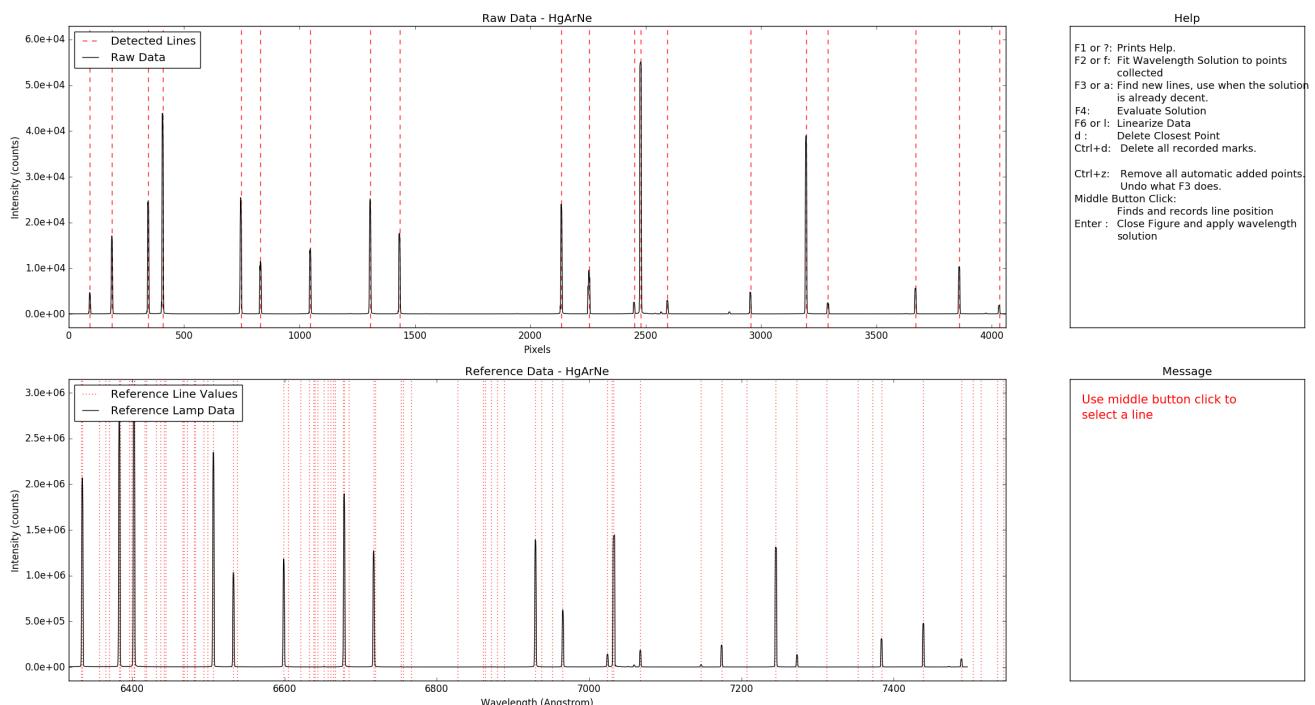


Running the Pipeline

2. The way you interact is by using your mouse and keyboard. The following table describes all the available keyboard commands with their respective keys and/or buttons, and a summary of their functionality.

Main Key	Alternative Key	Mouse Equivalent	Description
?	F1	None	Print Help message on the terminal
f	F2	None	Fit a function to collected points
a	F3	None	Find lines automatically
None	F4	None	Evaluate wavelength solution
d	F5	None	Remove point closest to the mouse pointer
l	F6	None	Linearize and smooth spectrum
m	None	Middle Button	Register the line closest to the mouse pointer
ctrl+z	None	None	Deletes ALL automatically added points
ctrl+d	None	None	Deletes all recorded points
ctrl+q	None	None	Ends the program
Enter	None	None	Accepts wavelength solution and closes the figure

It is advisable to practice a little bit to get familiar with this combination of keys and functions. Since we are still in the development stage, if you feel that the use of a key creates problems for you let us know.


In the following image the mouse pointer was placed very close to the center of the line and then with a middle button click the line is marked with a red circle. The software will calculate a centroid (vertical dashed blue line) and then will try to find the closest reference value (vertical red line). Keep in mind that the reference data (red dotted lines) are not the *reference lamp* lines themselves but values obtained from the NIST site indicated above.

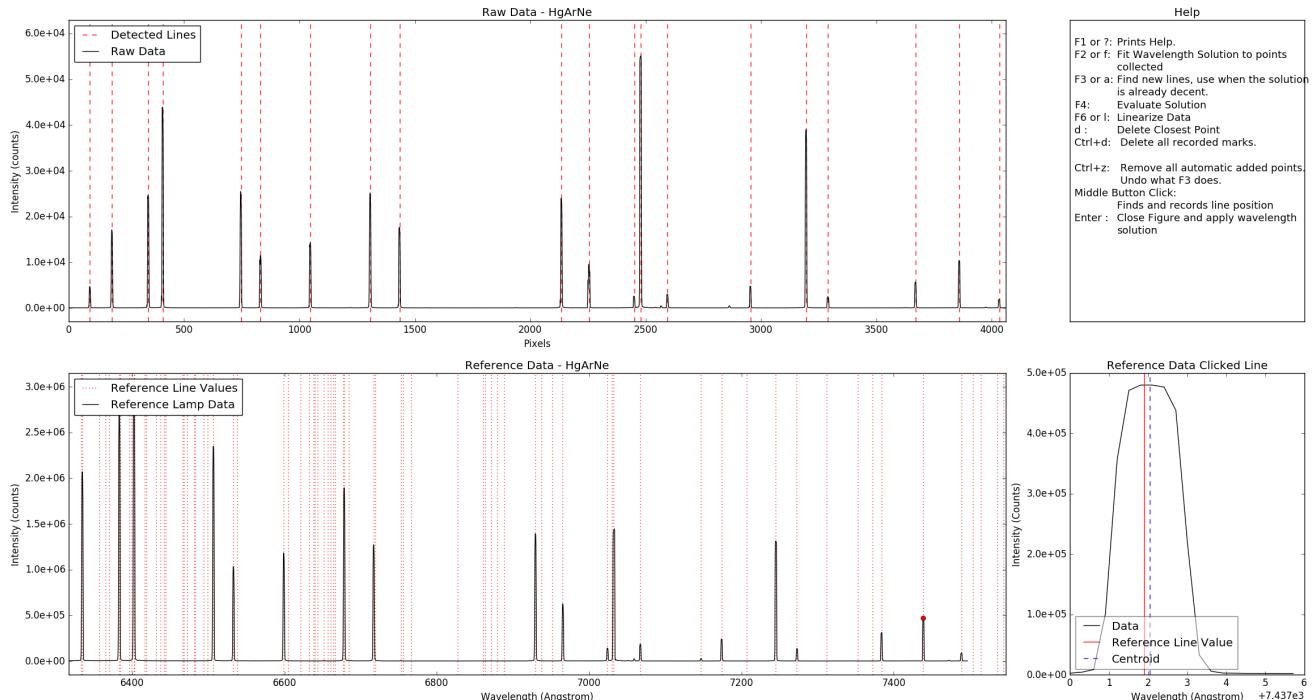


Then you have to find the corresponding line in the opposite plot. There is no preferred plot to start as long as you are consistent with the order in which you mark the lines.

3. In case you miss-identified a line, like in the image below, you can place the mouse pointer above the corresponding red circle and press **d** to delete it. It will search for the closest mark along the horizontal axis and remove it. If there is a counterpart in the opposite plot it will delete it too.

Running the Pipeline

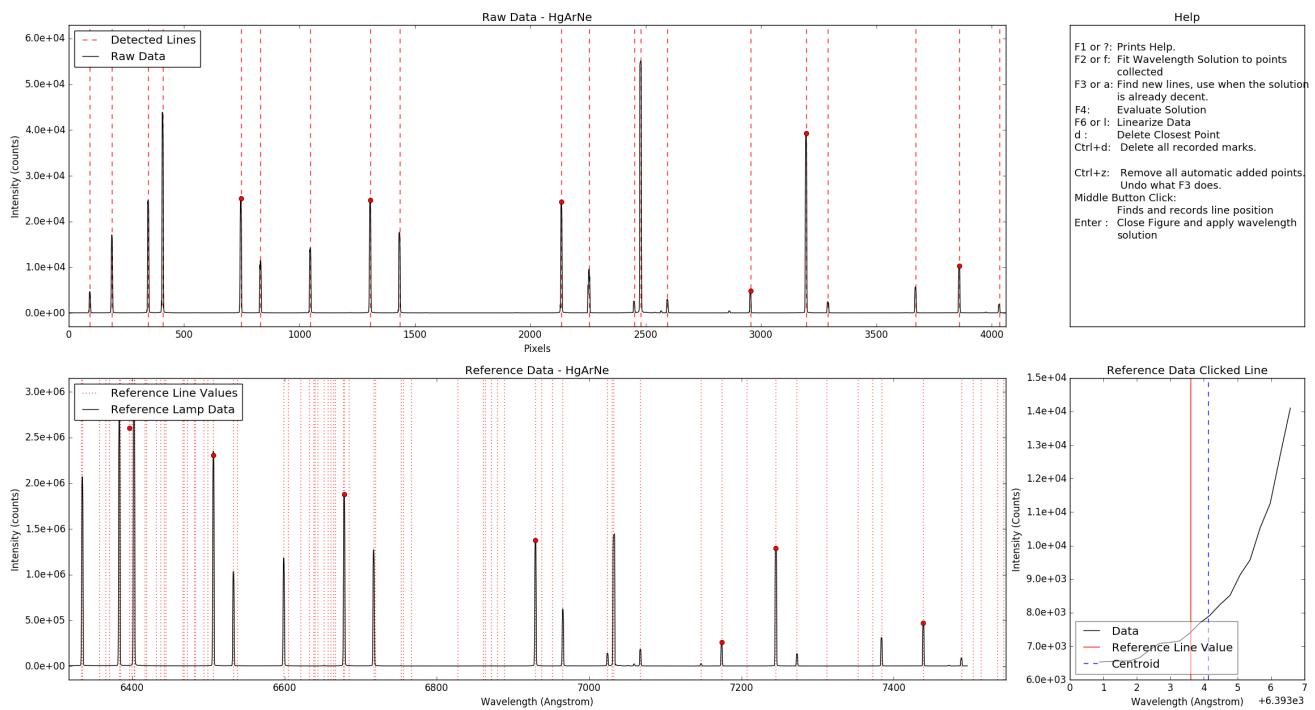

4. Once you matched a good number of lines, the minimum required by the fitting routines are four, you can either press F2 or **f** to make a fit of the pixels and angstrom values collected. Now the *Lower Right panel* will show the scatter plot of the fit. It is important to note here that the points in this plot do not represent the points you marked but rather each of the lines detected in your extracted 1-D comparison lamp spectrum (red dashed lines in the *Upper main plot*) It does one iteration of a 2-sigma clipping to reject outliers, and then it uses those values to calculate the Root Mean Square Error. In the example we show here the RMS is a bit high, but we will fix that below.

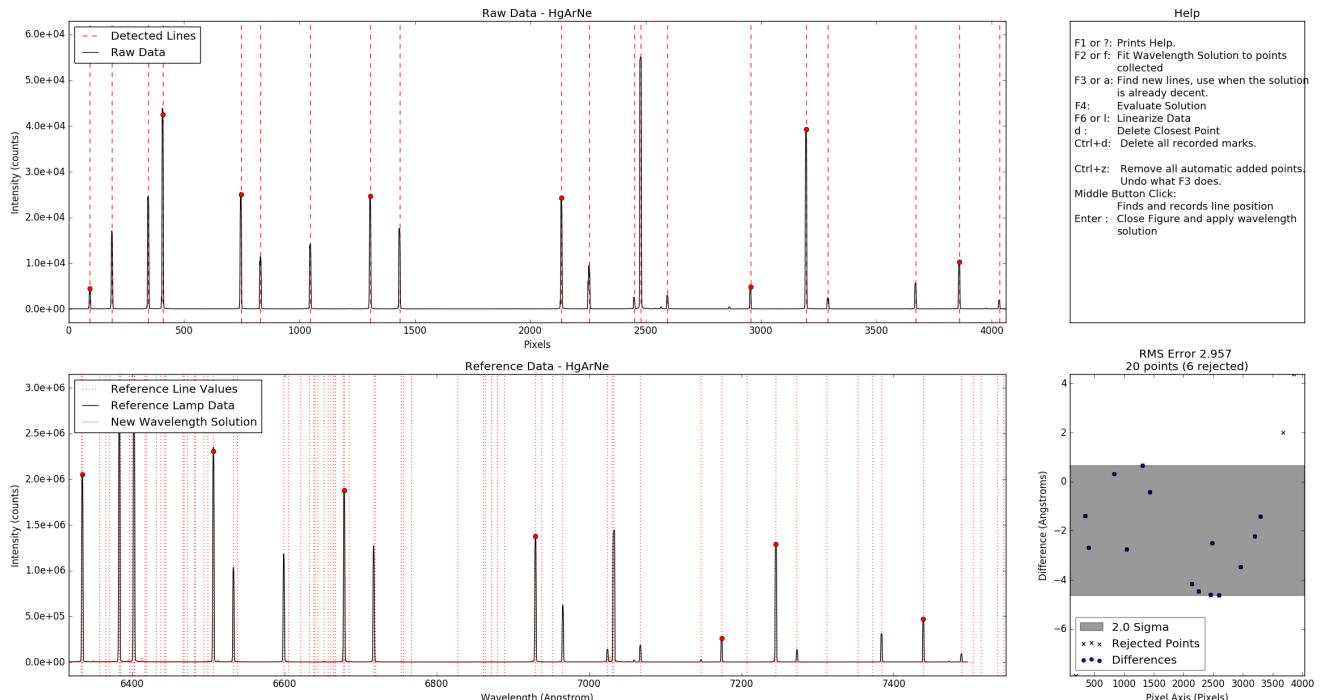


5. If you see that your current solution is decent you can press F3 or **a** and the pipeline will try to find more points automatically. The new matches found


by the software are shown as red dots at the base of each line in the two main *Upper* and *Lower* panels.

The automatic finding routine is not perfect, and indeed it depends on the preliminary wavelength solution. It uses the detected lines in the uncalibrated 1-D lamp, applies the preliminary solution and tries to find a match in the reference line values. In most cases it improves the solution, but not always so keep that in mind. In this case the RMS error is reduced by almost half, which is good, but if you look closely you can see the mismatches; also the *Bottom right* panel will show you this.

Running the Pipeline

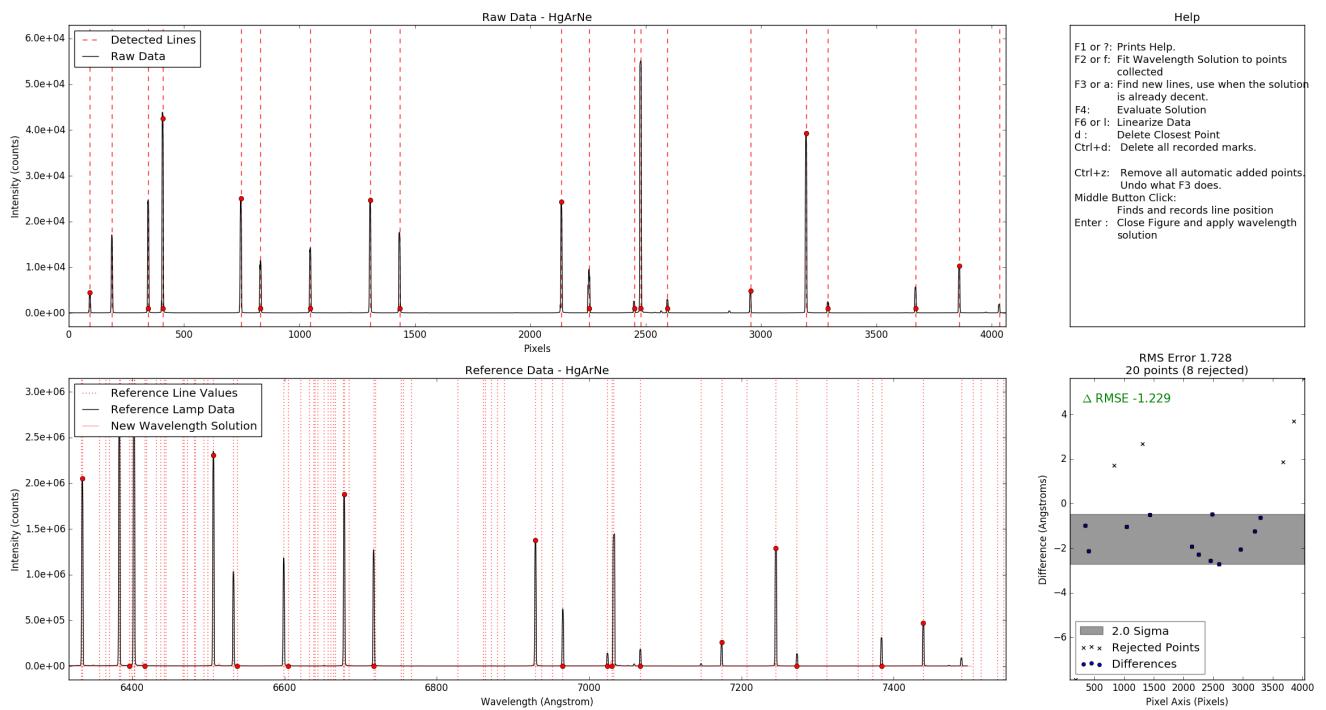


6. Now you we show an example of when the match is apparently good but in fact it's not. Here you need to zoom-in to see that there is an offset between the line centers and the matched laboratory line values, as shown in the figure below. You may have to apply different zoom values to your lamp and the reference lamp to get the plot to look like show it here. The solution in this case is locate the offending line(s), and delete it(them) pressing **d**. Then do a new match by clicking with the middle mouse button on the lines in the laboratory/reference plot, and the respective line in your 1-D uncalibrated lamp.

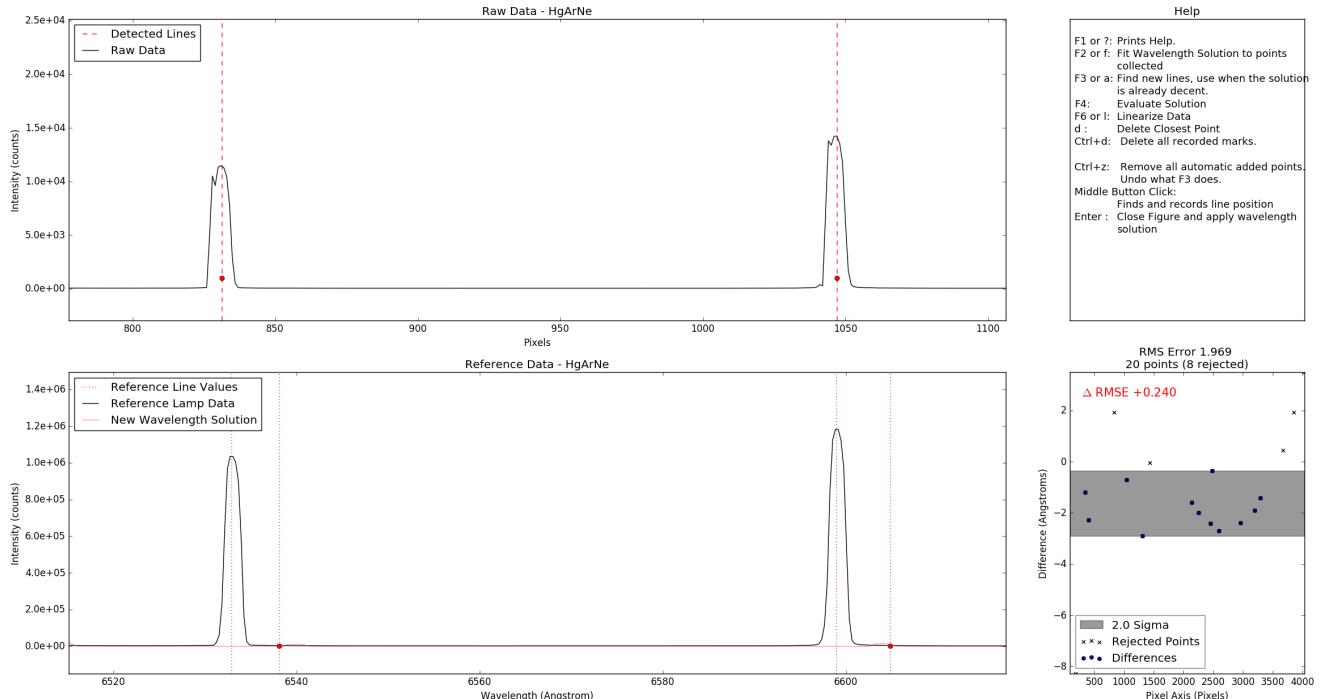


7. After you checked all the identifications and are happy with it, fit the solution again and you will obtain something like this:

Running the Pipeline


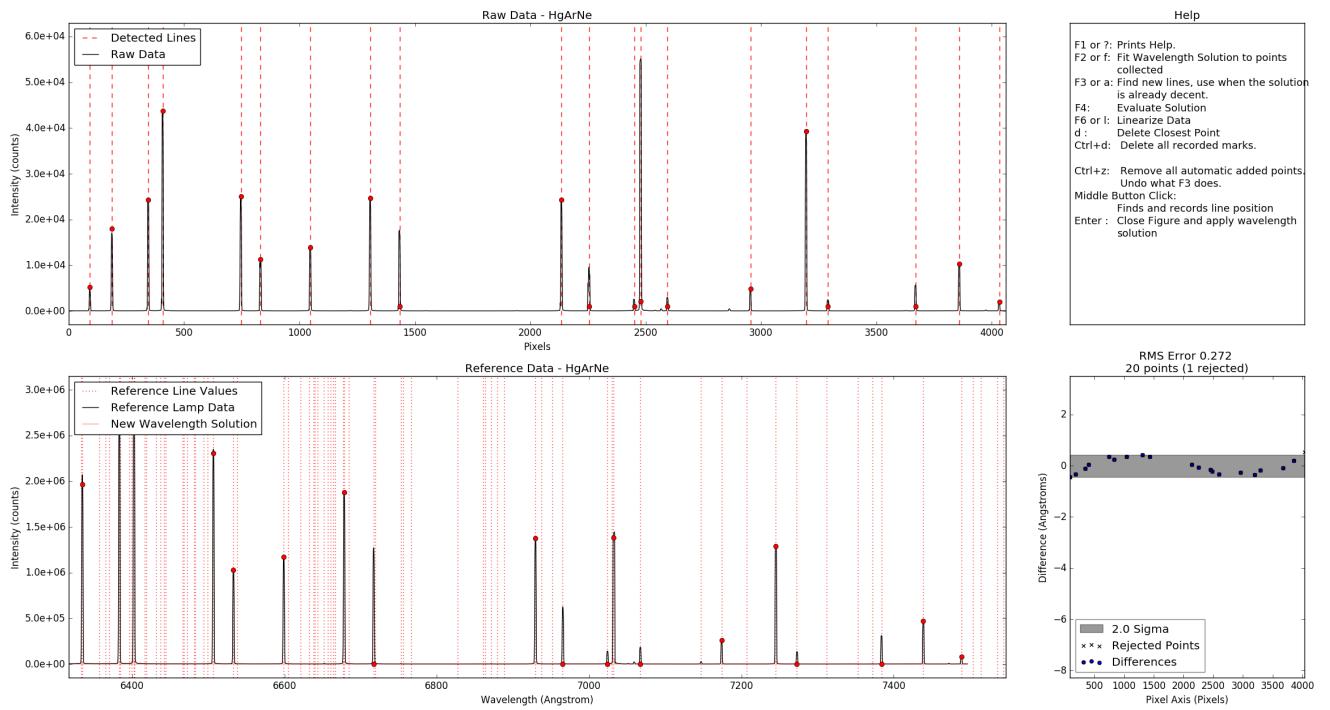


Here you can see that the *Bottom Right* panel shows the differences have a sinusoidal shape, which is also a sign that the solution can be improved. There are ways this can be implemented to refine the fit even further, but this is at present deferred to a later version.



8. Finally, a few samples of the spectra extracted by the pipeline.

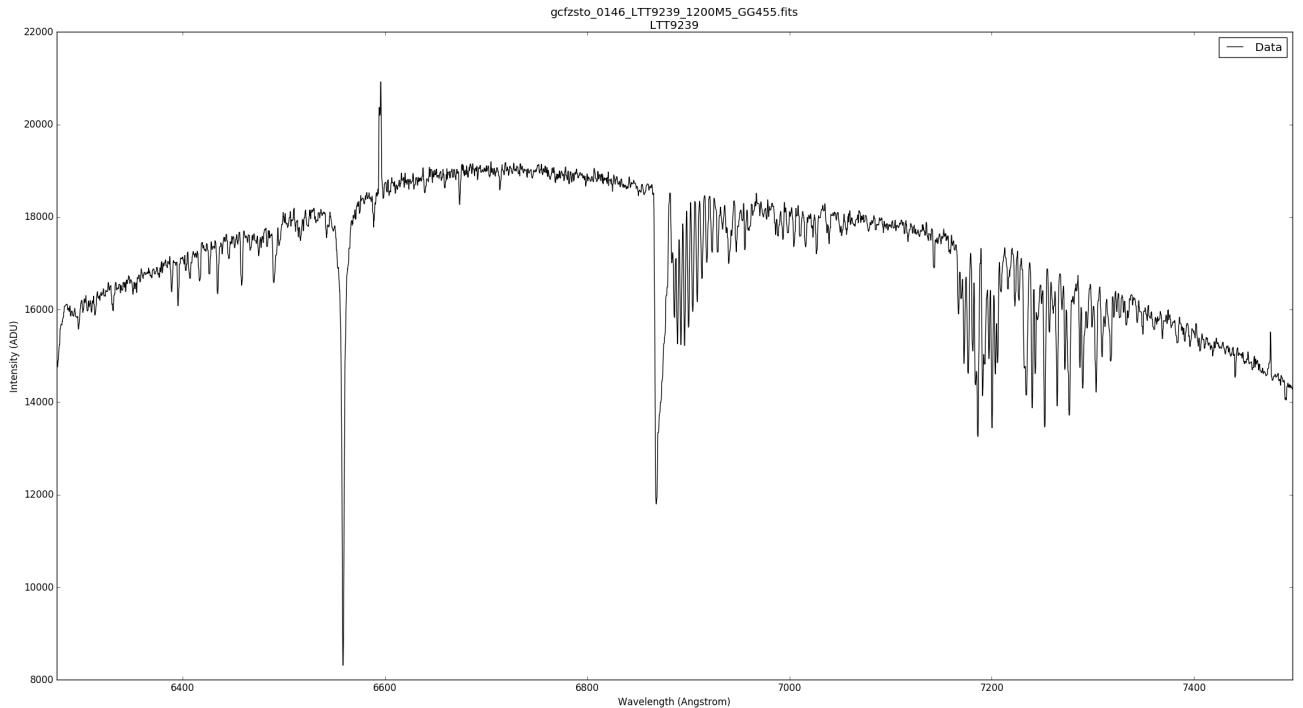


Running the Pipeline

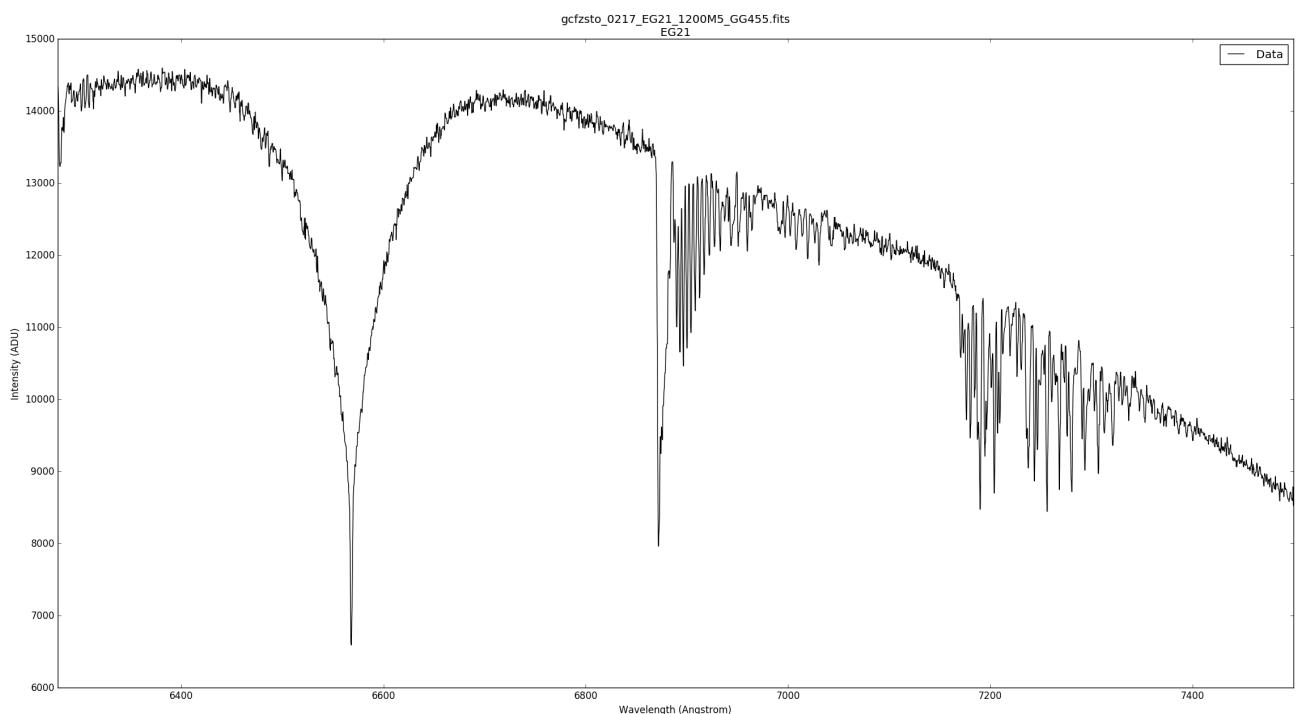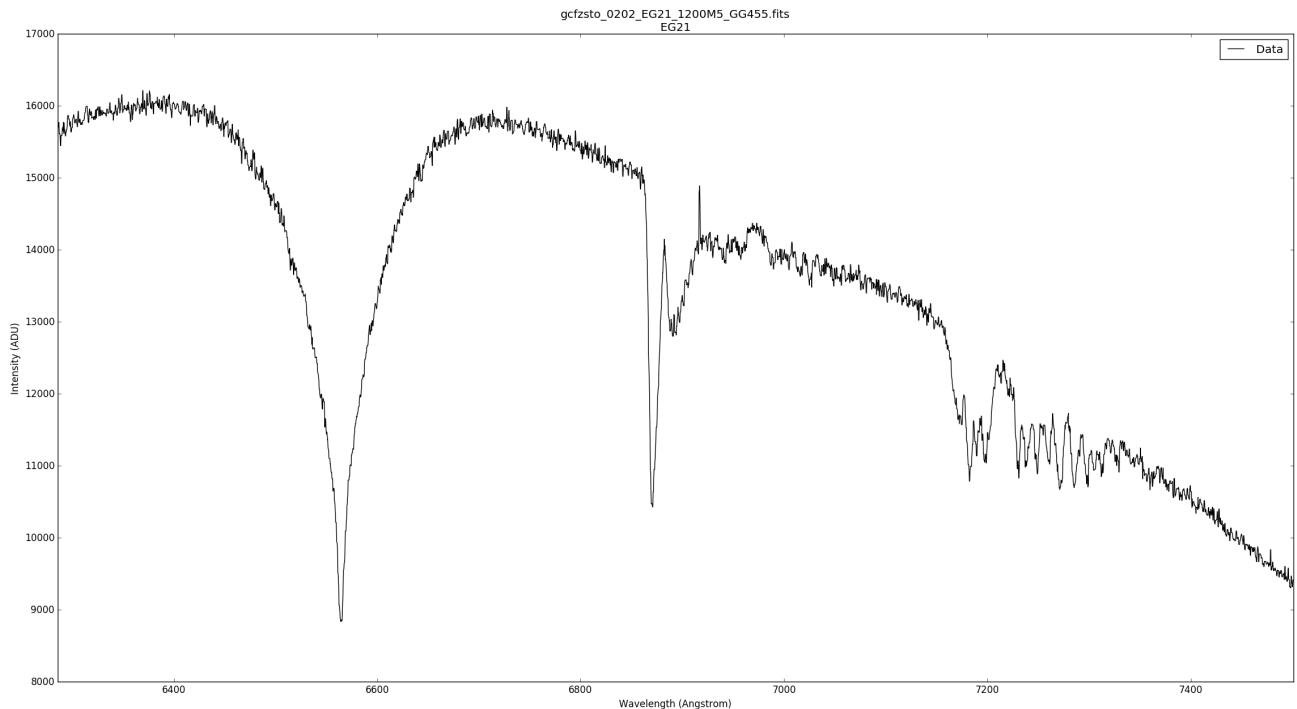

Running the Pipeline

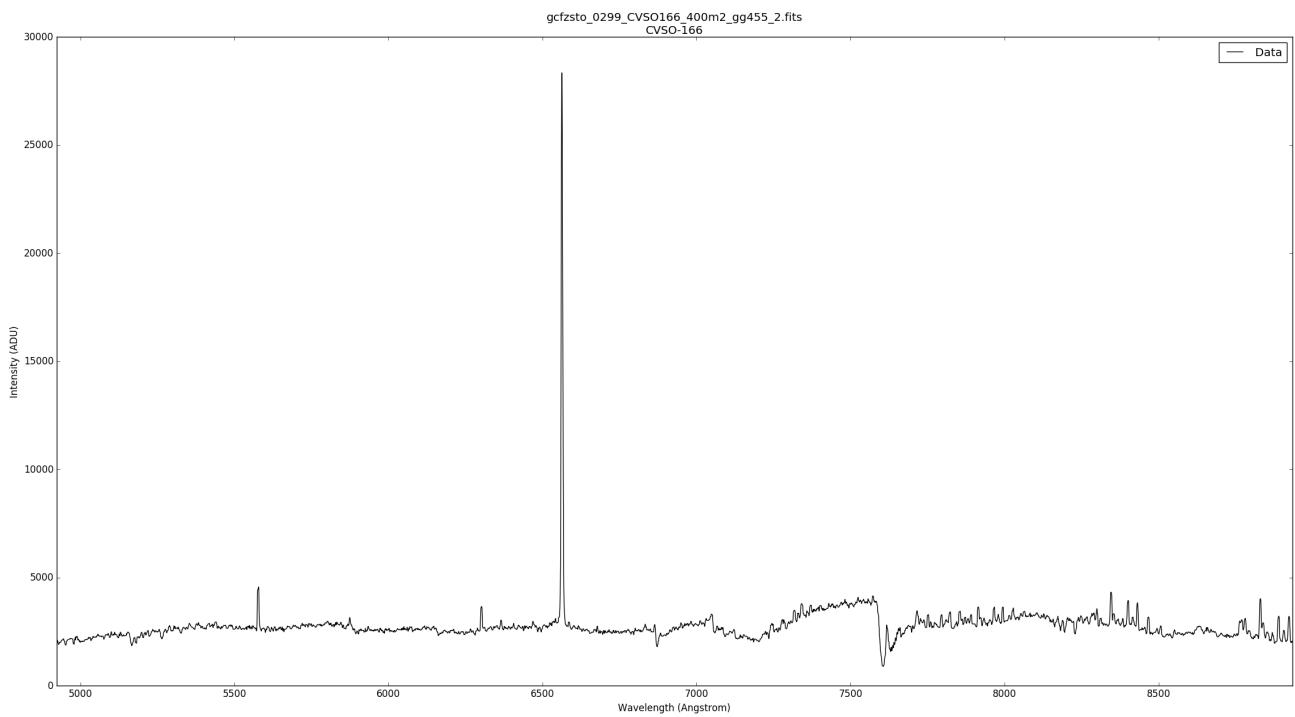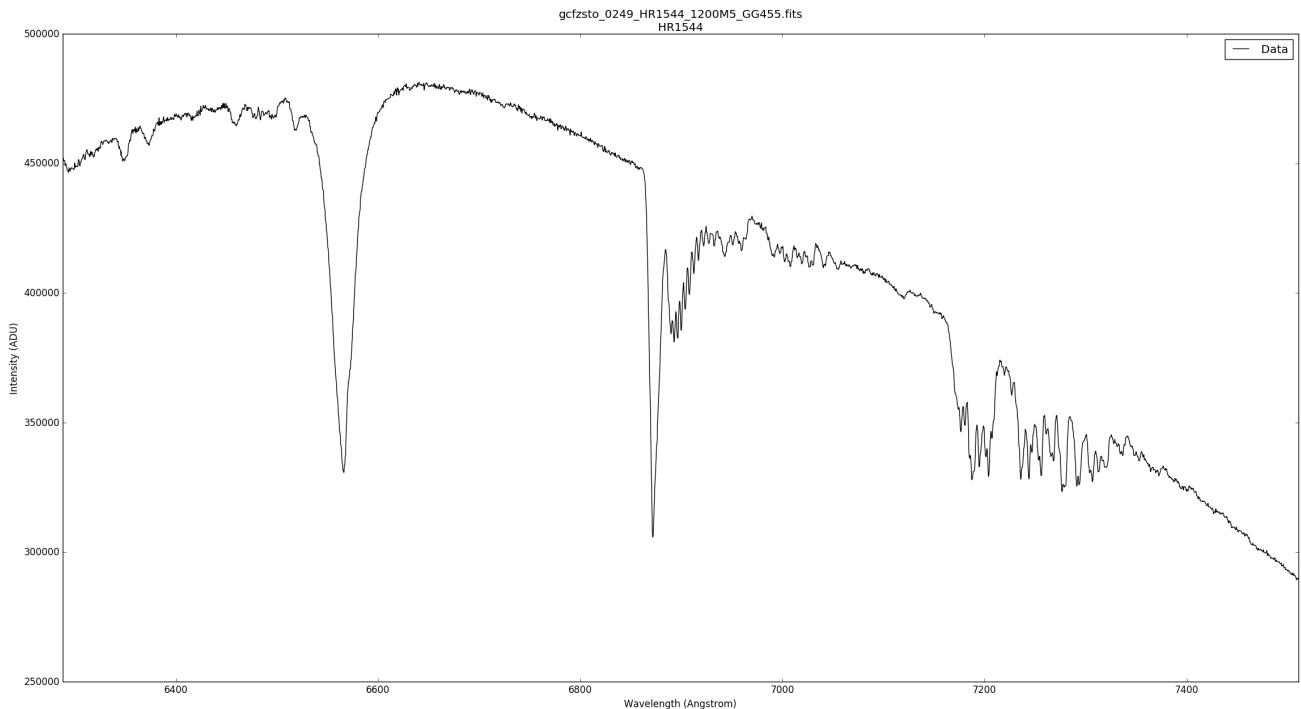

Troubleshooting

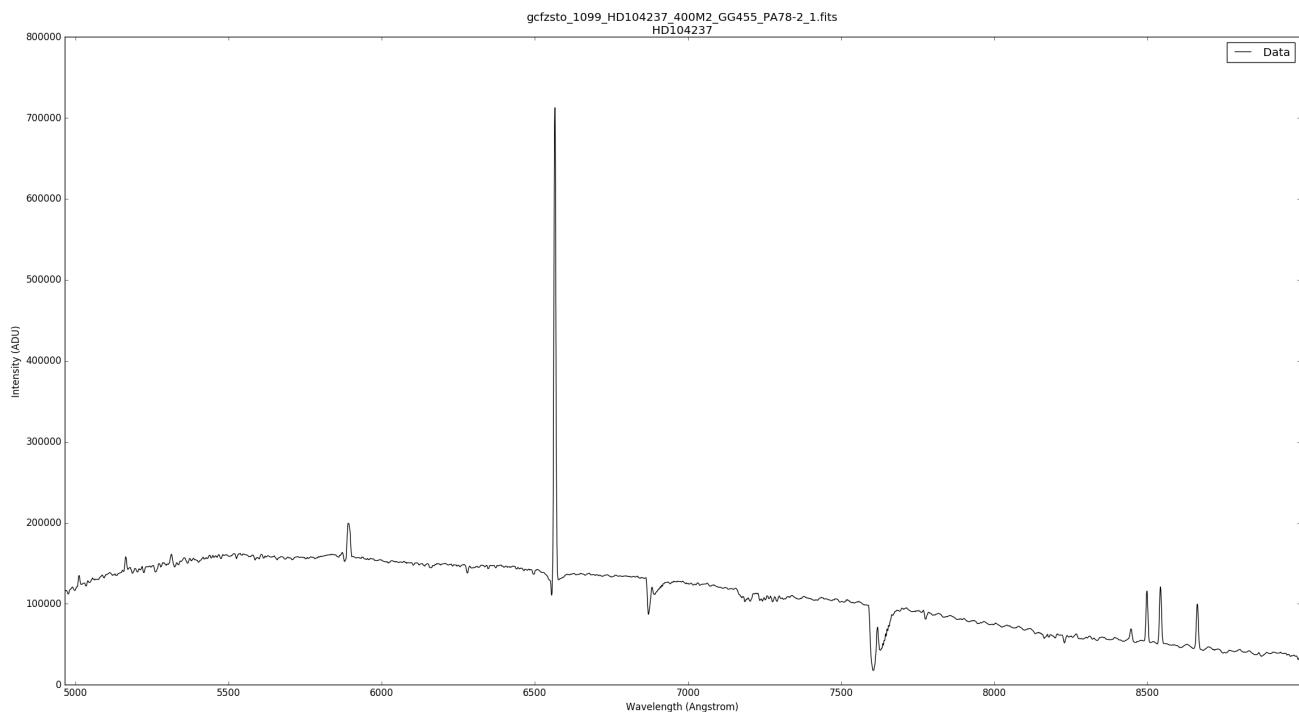

Troubleshooting

- The wavelength Solutions is way off: Check that the lamp was correctly registered in the header. Also check that the corresponding reference lamp exist. for instance is not the same to have HgArNe to HgAr
- Can't detect any objects: Check that the keyword `OBSSTYPE` is correct.
- The reference data plot, in interactive mode, doesn't show anything or only vertical dotted lines: The reference lamp doesn't exist for that configuration, since this is used only for visual reference sometimes it will display the same lamp but in other instrument configuration, this will not affect the quality of the solution.

Installation Instructions

Installation will slightly depend on the system but in general is simple and can be summarized in the following steps:

- Download the pipeline code
- Install prerequisites
- Install the pipeline

It is very important to note that this pipeline has been developed using *Python 2.7*, although we have done it in a way that would allow a smooth transition to *Python 3.5* we haven't tested it neither are we planning to do it in the near future.

Install Dependencies

There are two types of dependencies that have to be met. The system prerequisite installation depends on the platform itself so they will be detailed below in their respective subsection.

The python libraries are specified in the file `requirements.txt` and installing them is very easy. But it has to be done later on.

Ubuntu 16.04

Some other python-specific tools, if you already run python code most likely you already have them.

```
sudo apt-get install python-setuptools python-dev build-essential  
sudo easy_install pip
```

We have decided to use Qt4Agg backend since Qt seems to be the most multi platform compatible backend.

```
sudo apt-get install python-qt4
```

CentOS 7

Start by installing the EPEL repository

```
sudo yum -y install epel-release
```

Update the database with

```
sudo yum -y update
```

This takes a while...

Install pip

```
sudo yum -y install python-pip
```

Upgrade pip

```
sudo pip install --upgrade pip
```

Install python-devel

```
sudo yum install python-devel
```

Installing on MacOSX

Although I have successfully run the pipeline in two separated Mac OS X machines the installation process was a bit tricky and has not been fully tested. I plan to update this section in the next weeks.

Install Using Virtual Environments

Virtual Environment installation has not been fully tested.

Downloading the Goodman HTS Spectroscopic Pipeline

In order to get the code for the pipeline there are many options, our suggestion is to download the official release tar.gz file

<https://github.com/soar-telescope/goodman/blob/master/dist/goodman-1.0b1.tar.gz>

Installing the Pipeline

First of all install the python requirements. Your location must be the same as the file `requirements.txt` which should be your recently cloned repository

```
sudo pip install -r requirements.txt
```

Once this has succeeded proceed to install the pipeline using:

```
sudo python setup.py install --record files.txt
```

This will install the pipeline in your system and also will create a file `files.txt` that contains the list of files created at installation time and will be very helpful if you ever want to fully remove the pipeline.

Install DCR

In terms of cosmic ray rejection we shifted to a non-python package because the results were way better compared to LACosmic's implementation in astropy. LACosmic was not designed to work with spectroscopy though.

Visit this [Link](#) to download the code and find the instructions for compiling. I have added a few pre-compiled binaries and if you are lucky they will work right away. The available binaries are located in `goodman/dcr` and the options are:

- dcr.Ubuntu16.04
- dcr.Centos7
- dcr.MacOSSierra
- dcr.Solaris11

Choose whatever version fits your needs and rename it `dcr` and put it in a folder that at the same time is in your `$PATH` variable. If you don't know what that is follow the next section.

Install binary DCR

1. Open a terminal
2. In your home directory create a hidden directory `.bin` (Home directory should be the default when you open a new terminal window)

```
mkdir .bin
```

3. Move the binary of your choice and rename it `dcr`. If you compiled it most likely it's already called `dcr` so you can ignore this step.

```
mv dcr.Ubuntu16.04 ~/.bin/dcr
```

4. Add your `$HOME/.bin` directory to your `$PATH` variable. Open the file `.bashrc` and add the following line.

```
export PATH=$PATH:/home/myusername/.bin
```

Where `/home/myusername` is of course your home directory.

5. Close and reopen the terminal or load the `.bashrc` file.

```
source ~/.bashrc
```